

STANDARD PAPER SHEET FORMER

To manufacture paper samples in the Laboratory

According to: ISO 5269/2 - DIN 54.358 parte 1 - UNE 57.042/2 - NF Q 50.002.

In order to knowing the physical-mechanical and optic properties of the paper, cardboard and similar materials, it is necessary to manufacture them in the laboratory with a method quick and reliable that allows a sequence of quick realization of leaves under some similar conditions

The equipment is formed :

- A base in communication with the formation recipient
- A grill covered with a metallic mesh
- A cylinder of transparent plastic. This cylinder is adjusted to the circular grill
- An aspiration system in form of siphon.

The accessories have the strictly defined characteristics for the standard

- Stirrer (fig.1)
- Couch roll (fig.2)
- Tool for metallic mesh change (fig.3)




(Fig.1)



(Fig.2)



(Fig.3)

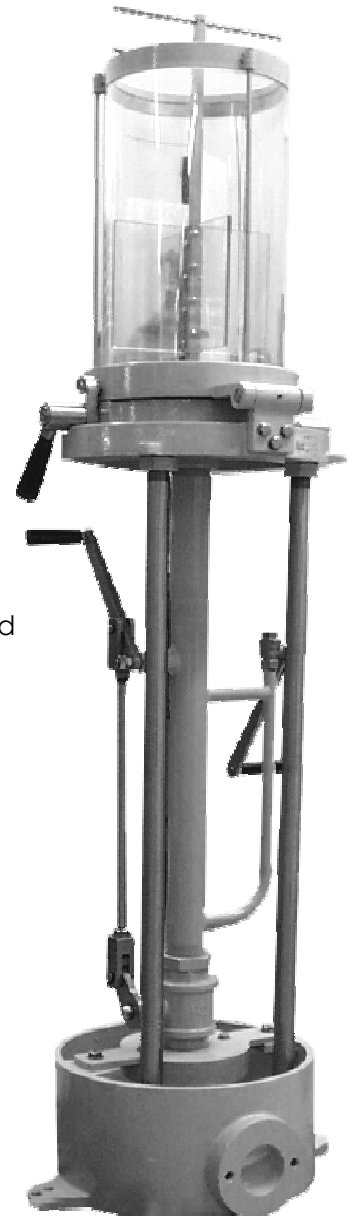

Height (m)	Width (m)	Depth (m)	Weigh (kg)
1.35	0.45	0.86	67