


WHITE WATER RECIRCULATION SYSTEM

F/ Standard Sheet Formers

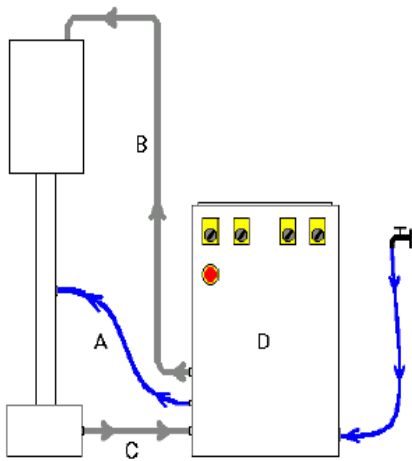

The equipment that is described below, when it works with a sheet former allows to reuse the used white waters when making standard samples of paper sheet.

The basic functions that can be controlled from the system are:

- Water inlet (A)**
- White Water Contribution (B)**
- White waters collection (C)**
- White Water Recirculation (D)**

Once the equipment is installed with the former, the appearance is the following:

To the left, the column where the paper sheet will be made, and to the right the white waters tank that includes the drive control, which represent graphically the functions that executes:

